



Just think of it as mission control.
Emphasis on control.




SciCan's variable speed E-STATIS™ electric motor system delivers near 'instant-on' control with speeds from 100 to 40,000 RPM and one of the shortest motors in the industry. The remarkable touch screen control (detachable for convenient placement) offers a range of presets and customization options. And, leave it to the Infection Control people to develop a motor that is fully autoclavable!

For information on the E-STATIS please visit www.scican.com.

Your Infection Control Specialist™

SciCanDental

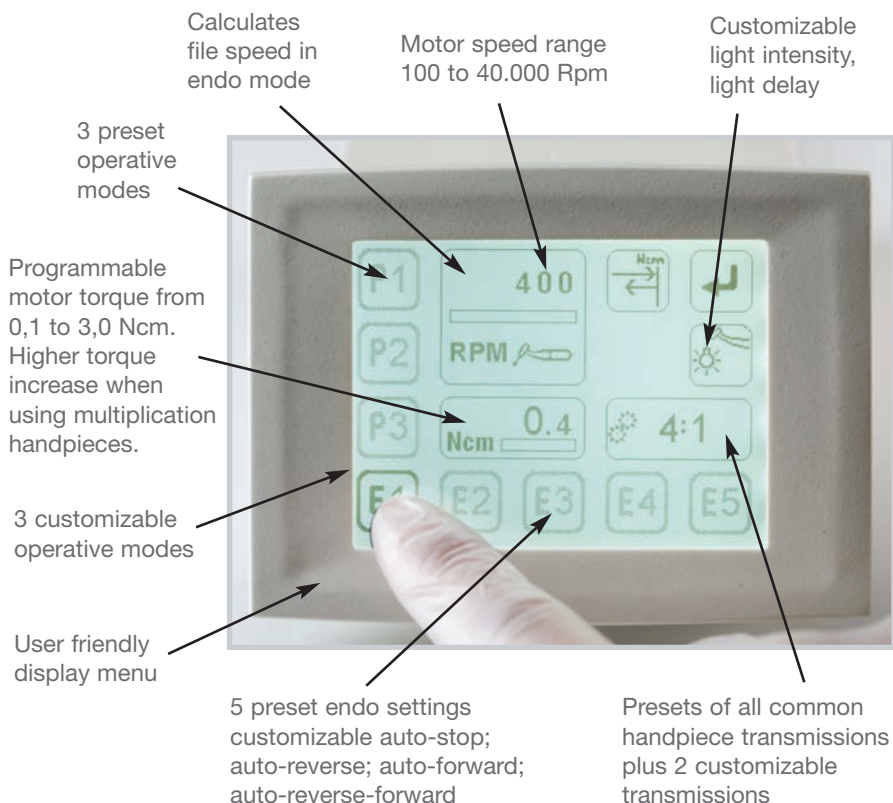


SPECIFICATIONS:

Motor speed range:	100 - 40.000 min-1
Motor power:	max 47 W
Motor voltage:	max. 22 V AC
Output torque:	max. 3,0 Ncm
Operation voltage of the bulb:	3.2V DC
Power of the bulb:	2.5W
Air outlet of the coupling (cooling air)	7 - 10 NL/min
Spray air pressure:	1 - 2.5 bar / 14.5 - 36.2 psi
Spray water pressure:	0.8 - 2.0 bar / 11.6 - 29.0 psi
Weight	99g
Dimensions control unit:	4.33 x 2.99 x 6.89 inch (110 x 76 x 175 mm)
Dimensions motor:	1.65 x 0.87 inch (42 x 22 mm)
Rotation:	clockwise / counter clockwise
Connection of attachments	ISO 3964

PRODUCT CODE:

S21100	E-STATIS Electric Motor System
S21103	E-STATIS Display extension kit

Manufactured by:
SciCan GmbH
 Wangener Strasse 78
 88299 Leutkirch im Allgäu
 GERMANY
 Tel.: +49 (0)7561 98343 - 0
 Fax: +49 (0)7561 98343 - 699

 0297

SciCan Inc.
 701 Technology Drive
 Canonsburg, PA 15317
 USA
 Phone: 724 820 1600
 Toll free: 1-800 572 1211
 Fax: 724 820 1479

SciCan Ltd.
 1440 Don Mills Road
 Toronto ON M3B 3P9
 CANADA
 Phone: 416 445 1600
 Toll free: 1-800-667-7733
 Fax: 416 445 2727

Three year motor warranty.
 Fully autoclavable motor.
 Brushless technology – greatly reduces service requirements.
 Short and lightweight motor – weighs only 99g.



Compact design for easy installation



Accurate calculation of bur & file speed display



Easy retrofitting of pneumatic dental chairs.
 Power plug can be taken off the cable completely for easy integration!