



**It's 99.9999% deadly.
Just not to you.**



Virucidal. Bactericidal. Tuberculocidal. Just not harmful to you or your patients. OPTIM® disinfecting wipes kill germs on surfaces fast – up to 10 times faster than other leading cleaners. OPTIM cleans supreme using a patented formulation based on Hydrogen Peroxide that has virtually no odor. Also, the solution readily biodegrades into water and oxygen after disinfection. So OPTIM is eco-friendly and people friendly. In fact, it's really only germs that aren't too fond of it. Take control, because the stakes are too high.™


Your Infection Control Specialist™

SciCanDental



OPT/IM 33TB WIPES – EXCEPTIONAL CLEANING & DISINFECTION CAPABILITIES

- One step disinfectant cleaner – fast & convenient
- AHP – patented technology designed to deliver better cleaning and faster disinfection results (active ingredient 0.5% hydrogen peroxide)
 - Hydrogen peroxide does not leave behind active residue which may lead to the development of resistant strains of microorganisms
 - OPT/IM Wipes remain extremely wet and do not dry out

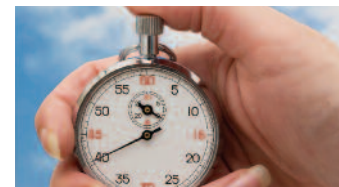

FAST CONTACT TIME:

- 1 minute bactericidal & virucidal, tuberculocidal
- 30 second sanitizer

SUPERIOR ENVIRONMENTAL AND SAFETY PROFILE

- No strong chemical odor
- Active ingredient breaks down into oxygen & water within 1 minute of disposal into sewage system
- Non-flammable
- Non-irritating to skin & eyes
- 0 HMIS hazard rating*, non-toxic, EPA Cat. 4, the lowest toxicity level on EPA toxicity category

Legend HMIS	
Severe	4
Serious	3
Moderate	2
Slight	1
Minimal	0

Health	0
Flammability	0
Physical Hazard	0

* Hazardous Materials Information Systems – a requirement on MSDS by OSHA

COMMON ACTIVE INGREDIENT & AHP

	PHENOLS	QUAT(S) ALCOHOLS	AHP (OPT/IM 33TB)
Toxicity:	Highly toxic (o-phenylphenol recognized to cause cancer in the state of California – Proposition 65)	Contain hazardous ingredients	EPA toxicity category 4 – no warnings required on label. Non-toxic (Oral LD > 5000mg/kg)
	Irritating to eyes & skin	May cause dry skin & eyes and skin irritation	Non-irritant to eyes & skin
	May contain hormone disrupting agents	Contribute to poor indoor air quality (VOC's)	No noxious fumes
Hazard Facts:	Combustible	Volatile, flammable and must be stored in well-ventilated area	Classified as non-hazardous and non-flammable
		Poor biodegradability	No waste treatment required, biodegradable
Contact Times:	Long contact times (e.g. 10 minutes)	Alcohol presence allows for relatively short contact times	Short contact times (1 minute virucidal & bactericidal, tuberculocidal)
Precautionary statements examples:	“DANGER” “severely irritating to eyes and skin, combustible” “hazard to humans and domestic animals”, corrosive, causes irreversible eye damage”	“CAUTION” “hazard to humans and domestic animals”, “causes eye irritation”, “probable mucosal damage may contraindicate the use of gastric lavage”, “Flammable”	Non-irritant to eyes & skin according to criteria available under OECD. No warnings except for the standard statement “Keep out of reach of children!”
Active Ingredient %	Up to 15%	QUAT-Alcohol: 10-25%, Alcohols: 60-80%	0.5% Hydrogen Peroxide

ORDERING INFORMATION FOR OPT/IM 33TB

Dry Wipes:

SCI33R18 Optim 33 TB Dry Wipes Refill 6 x 7”
(1 case: 18 x 160 sheets)

SCI33R18XL Optim 33 TB Large Dry Wipes Refill 10 x 10”
(1 case: 18 x 60 sheets)

Liquid:

SCI3312X1-ALOE Optim 33 TB with Aloe Scent, 1L (1 case: 12 x 1L)

SCI334X4 Optim 33 TB unscented, 4L (1 case: 4 x 4L)

SCI334X4-ALOE Optim 33 TB with Aloe Scent, 4L (1 case: 4 x 4L)

SciCan Ltd., 1440 Don Mills Road, Toronto, Ontario, M3B 3P9 • Phone: 416.445.1600 • Toll free: 1.800.667.7733 • Fax: 416.445.2727

SciCan Dental

Your Infection Control Specialist™